




Details zu HPLC-Applikation Nr. 116

Nucleoside auf CS-ODS-5 μ

Säule: 150 x 4,6 mm CS-ODS 100-5 μ

Art.-Nr.: 586 1310 15

HPLC-Bedingungen: Eluent: Wasser
Fluss: 1,3 mL/min
Temperatur: RT
Detektor: UV 254 nm

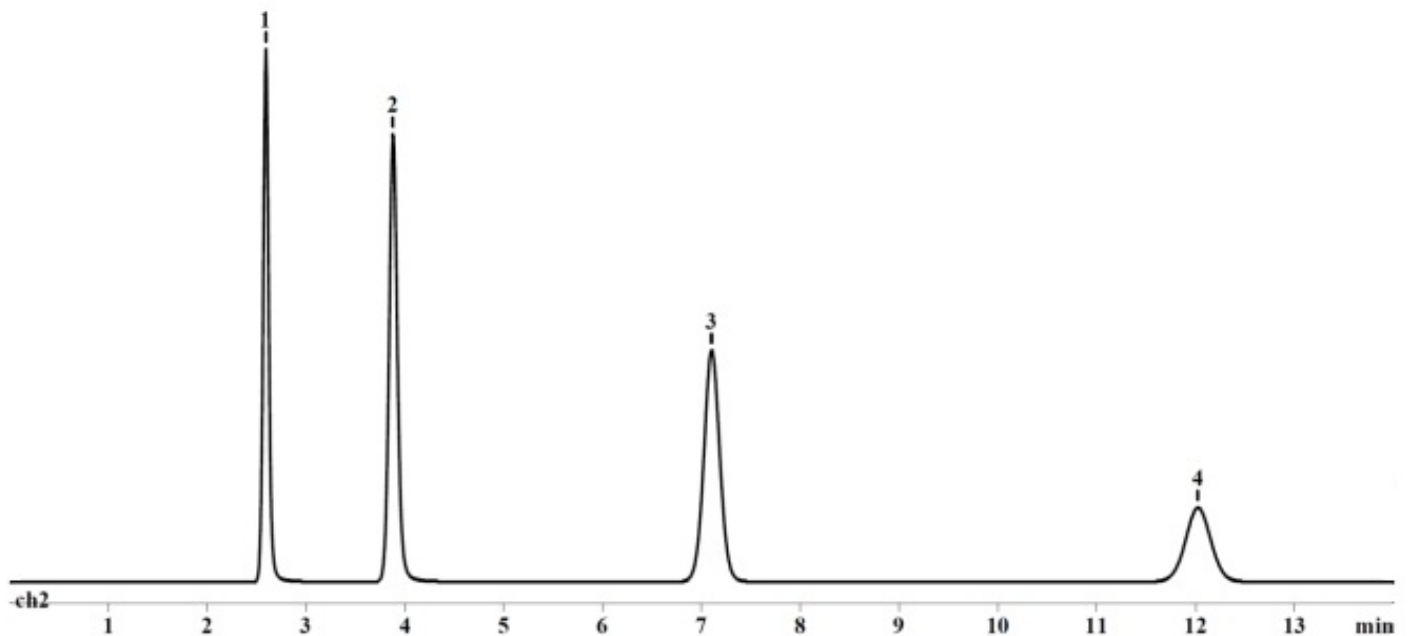

1.) Cytosin

2.) Uracil

3.) Cytidin

4.) Uridin

Peaks: Cytosin
Uracil
Cytidin
Uridin